

TARTAR & STAIN



Removes calculus, food stains, permanent cement, tartar and tobacco

SX4061

3.8litre

TEMPORARY CEMENT REMOVER



Especially formulated for most Zinc oxide/eugenol temporary cements

SX4062

3.8litre

PLASTER & STONE REMOVER



Removes plaster and stone from partials, flasks and articulators

SX4063

3.8litre

BRACON JEWELLERY CLEANER



General purpose non ammoniated. Concentrated 1:7 (makes 31.5 litres of solution).

SL1002

4.5litre

BRACON CLINICAL DETERGENT

General purpose non-ammoniated. Concentrated 1:20. (makes 100 litres of solution).

SL1001

5litre



L&R TEMPORARY CEMENT REMOVER

Ready to use.

SL1005

3.8litre

L&R PLASTER + STONE REMOVER SOLUTION

Ideal to remove plaster and stone from dentures, partials, flasks, articulators and general instruments in an **ultrasonic** unit. Ready to use.

SL1012

3.8litre

L&R GENERAL PURPOSE CLEANER

Non-ammoniated. Removes blood, tissue, debris and compounds from instruments, burs, general equipment, orthodontic appliances, gold and porcelain. Concentrated 1:10 (makes 38 litres of solution).

SL1011

3.8litre

L & R TARTAR + STAIN REMOVER LIQUID

Cleans tartar, light stains and cement from dentures, partials and bridges. Ready to use.

SL1004

3.8litre

CLEANER	CONTAMINANT	APPROX CLEANING TIME	SPECIAL NOTES
Tartar, Light Stain and Permanent Cement Remover SX 4061 SL1004 SL1021	Tartar & Light Stain Permanent Cement Solder flux Oxides and Investments Rust Hard water scale	5-10 mins 1-10 mins 1-5 mins 5-10 mins 1-2 mins 1-2 mins	This is a non-fuming acid. When used on dentures, rinse, brush lightly and change after each use. If instruments remain more than 2 minutes, they will discolour (may be removed with a buffing wheel). Rinse in water and neutralise in non-ammoniated. Use in glass beaker only. USE FULL STRENGTH.
Temporary Cement SX4062 SL1005	Most Zinc Oxide/Remover Eugenol types	5-10 mins	Some temporary cements are only soluble in Tartar, Light Stain and Permanent Cement Remover. Rinse thoroughly with water. Change after each use. USE FULL STRENGTH.
Plaster and Stone Remover SX4063 SL1012 SL1025 SK7115 SC7019	Plaster & Stone from dentures, partials, flasks and articulators	5-10 mins	Rinse thoroughly with water. USE FULL STRENGTH.
Jewellery Cleaner, Non-Ammoniated SL1002	Buffing compounds, rouges, tripoli, oxides Dried blood and debris	1-5 mins	Dilute: 1 part cleaner solution to 7 parts water.
Clinical Detergent SL1001	Non-Ammoniated Bracon Clinical	1-5 mins	Dilute: 1 part cleaner solution to 20 parts water.
Non-Ammoniated General Purpose SX4060 SL1011 SL1022 SK7110	Cleans diamonds, burs, slabs, gold, instruments, porcelain etc. Rinse thoroughly with water.	1-5 mins	Dilute: 1 part solution to 10 parts water. Use in a glass beaker only. Rinse thoroughly with water.

SPECIAL NOTES CAUTION:

- Avoid skin contact with chemicals
- Do not use on natural teeth
- Use proper accessories for handling during cleaning procedure
- Rinse after every cleaning procedure to prevent staining
- Do not use on aluminium
- Not for internal use *N.B. For Ultrasonic Baths (SEE SECTION 1 PAGE 1.33)*

VACUUM FORMING BLANKS

BRACON COLOURED MOUTHGUARD BLANKS

These should only be treated as an approximate guide. Please phone for a colour guide.



For mouthguards. 4 mm thick please specify when ordering:
125 x125 Square or 120 mm ø Round

BRACON CUSTOM COLOURED




For mouthguards. 4 mm thick

SQUARE 125x125 mm

AN1001S	St George	each
AN1002S	St Andrew	each
AN1003S	3D	each
AN1004S	Union jack	each
AN1005S	Rainbow	each
AN1006S	Camouflage	each
AN1007S	Chequered flag	each

ROUND 120 mm ø

AN1001R	St George	each
AN1002R	St Andrew	each
AN1003R	3D	each
AN1004R	Union jack	each
AN1005R	Rainbow	each
AN1006R	Camouflage	each
AN1007R	Chequered flag	each

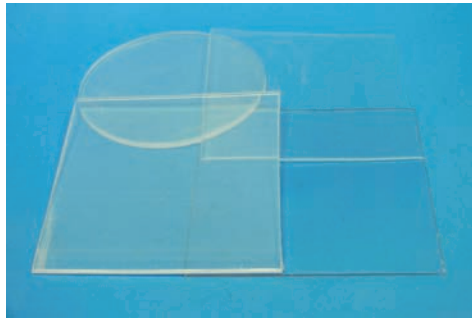

BRACON CLEAR SOFT BLANKS

For mouthguards. and bleaching trays

SQUARE 125 x 125 mm

SW1008	5mm	25
SW1009	5mm	5
SW1010	4mm	25
SW1011	4mm	5
SW1015	3mm	25
SW1015A	3mm	5
SW1020	2mm	25
SW1025	2mm	5
SW1018	1.5mm	25
SW1019	1.5mm	5
SW1016	1mm	25
SW1017	1mm	5

ROUND 120 mm ø

SW1010R	4mm	25
SW1011R	4mm	5
SW1015R	3mm	25
SW1015AR	3mm	5
SW1020R	2mm	25
SW1025R	2mm	5
SW1018R	1.5mm	25
SW1019R	1.5mm	5
SW1016R	1mm	25
SW1017R	1mm	5

BRACON CLEAR HARD BLANKS

Clear thin for splints and Clear or White thick for trays

SQUARE 125 x 125 mm

SW1033	0.5mm clear	10
SW1035	1mm clear	10
SW1037	1.5mm clear	10
SW1039	2mm clear	25
SW1040	2mm clear	5
SW1041	3mm clear	25
SW1042	3mm clear	5
SW1030	3mm white	25
SW1031	3mm white	5

ROUND 120 mm ø

SW1033R	0.5mm	10
SW1035R	1mm	10
SW1037R	1.5mm	10



LAMINATE SOFT/HARD BLANKS

Double layer, hard and soft for occlusal and grinder splints or thicker layers for functional positioners.

KOMBIPLAST



ROUND 120 mm ø

SP2103	1.8mm	10
SP2103A	3mm	10

ERKOLOC #591230



ROUND 120 mm ø

SN1010	3mm	10
---------------	-----	----



FORMING UNITS

DUNAFORM V



For vacuum forming base plates, bleaching trays, custom trays, temporary splints, mouth guards, bruxism splints, copings and more! Dunaform V saves time and money, is fast and easy to use.

- Extremely simple to handle
- Works independently of air pressure
- Both round and square foils with the dimensions \varnothing 120 and 125 mm and 125 x 125 mm can be processed
- For foil thicknesses of 0.1 to 4.0mm
- Heating from above
- Economically priced. One-year guarantee assured quality and value

Mains Supply	230V / 50Hz
Output	550w
Height	320mm
Width	260mm
Depth	260mm
Weight	7kg

ED1002 each

PYREX BEAKER



600ml capacity

SL1008 each

DUNAFORM PRESSURE



A pressure moulding device for all applications of the dental thermoforming technique in the laboratory and in the practice.

- The thermic-controlled IR-radiator gets to working temperature within a short time
- Heating and cooling times are programmable
- The heating of the sheets takes place on the model-turned-side
- Optical and acoustic warnings
- Large display-indication
- Simple operation
- Low air consumption levels

Mains Supply	230V / 50Hz
Output	850w
Height	185mm
Width	350mm
Depth	295mm

ED1001 each

Bracon is a proud supplier of:

DETAX

transactor

MESTRA®

BLUE FREEZE COLD SPRAY




Instantly cools the formed material to prevent pulling away from the model

SC3001 200ml

WORK TRAYS

MESTRA MODEL TRAY SYSTEM

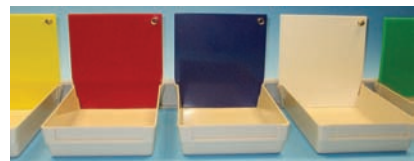

- Simple and elegant way of organising your work
- Provision for labels on tray front
- Colour coded back panel with clip for worksheet
- 5 colours – Yellow, Blue, Red, Green and White – one for each day of the week, or as you like

SM0100050	Yellow
SM0100052	Blue
SM0100054	Red
SM0100056	Green
SM0100058	White
SM0100060	Mounting rail